RECOVERY.M8

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
USER MANUAL

**Tak fordi du handlede hos os!
Thank you for shopping with us!**

Pumpeenhed

Display

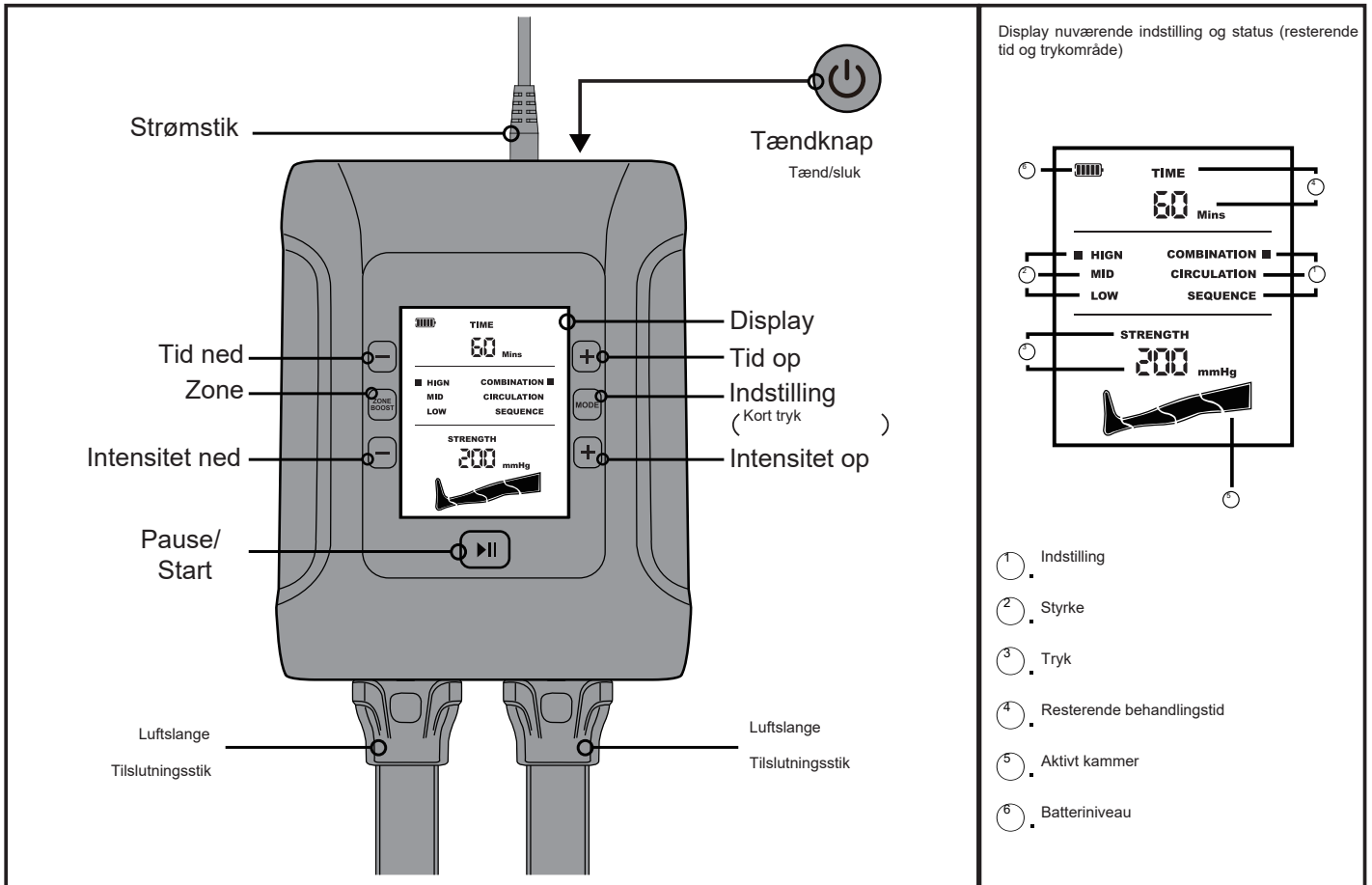Rengøring & Vedligeholdelse

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Warranty


The product comes with a full 1 year warranty. However, this can be waived if:

1. The product has been damaged as a result of irresponsible use or storage.
2. The product has been disassembled.
3. The instructions for use and cleaning in this manual have not been followed.
4. Lack of proof of purchase of the product (invoice).

Main unit

Display

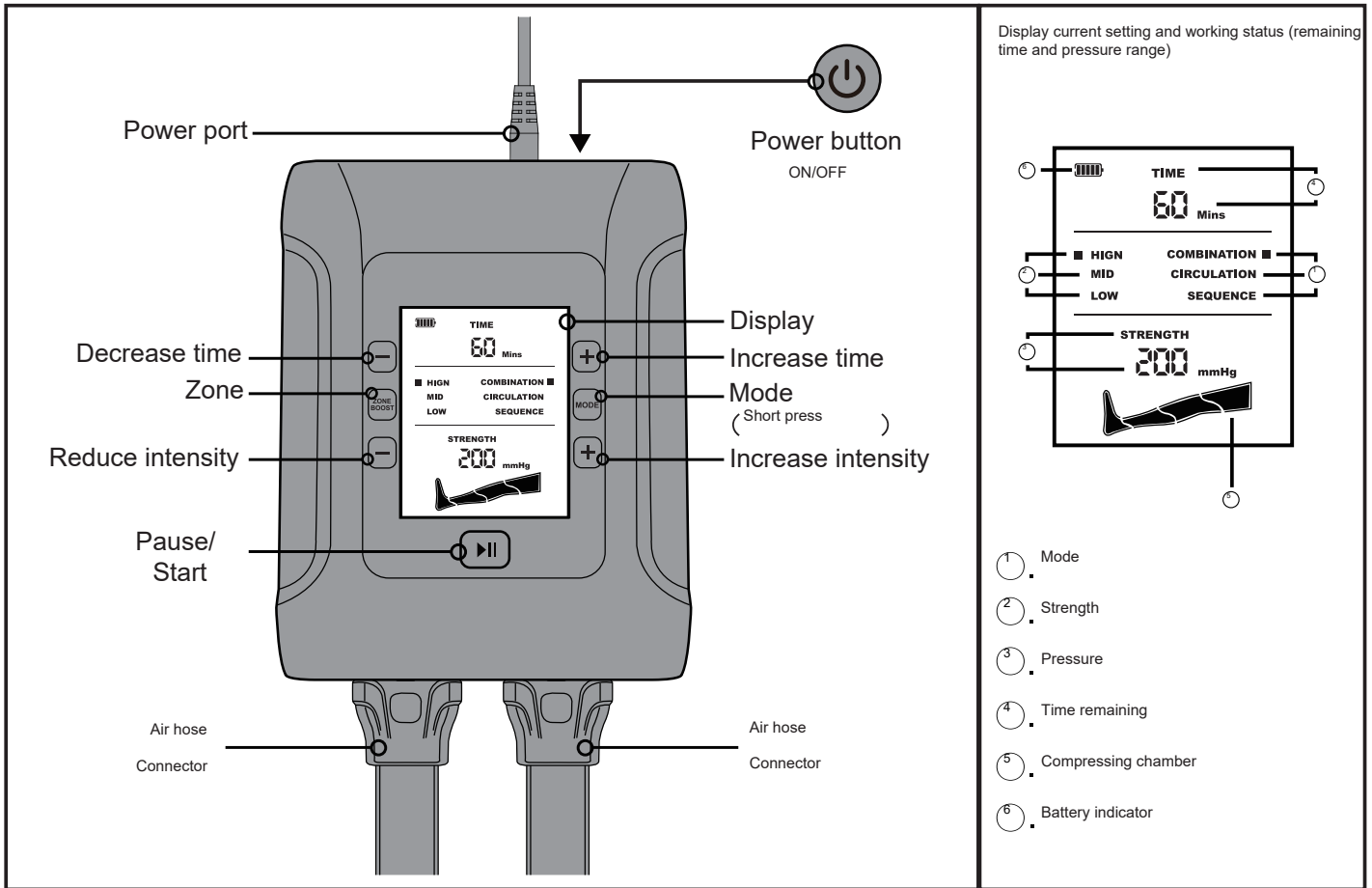Cleaning & Maintenance

- Switch off the pump before cleaning.
- Do not wash the pump or boots with water.
- Use a damp cloth to clean the pump and boots.
- Store clean and dry.
- Protect the pump from direct sunlight and other extreme heat sources.

Warning

- Do not unzip the boots while the pump is operating.
- Do not expose the pump to any liquid.
- Make sure the pump is on a flat surface.
- Do not adjust the boots while the pump is operating.
- The pump must not be repaired, maintained or cleaned when it is switched on.

This product must not be used in the following situations (not limited to), without the consent of your doctor. Acute inflammatory skin diseases, arrhythmia, erysipelas, deep vein thrombosis (DVT), pulmonary edema, uncontrolled hypertension, patients with pacemakers.



Quick start guide

Read the manual carefully before use. Follow the instructions to avoid malfunctioning of the product and to get the best experience.


1. Insert the sleeve hose connectors firmly into the hose.

* It is possible to use one sleeve at a time. If you want to do this, simply take the silicone plug on the back of the pump unit and insert it into the desired outlet.


2. Connect the hose to the control unit, you will hear an audible 'click'.

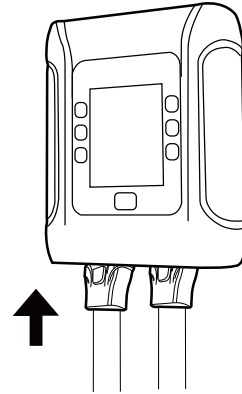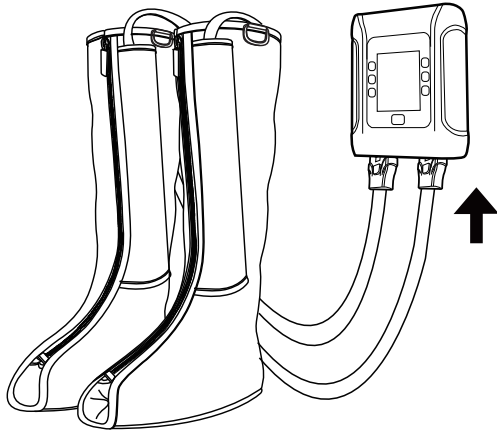3. Connect the adapter into the outlet and the control unit, the screen will show the charging battery indicator.

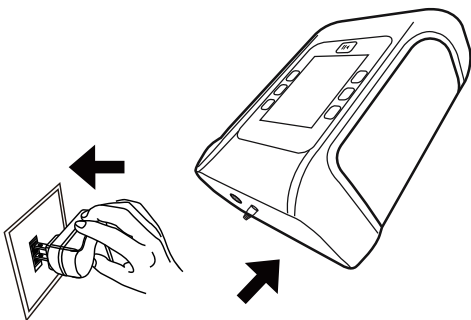4. Put on the sleeves correctly. Then set the desired program & strength and start the treatment.

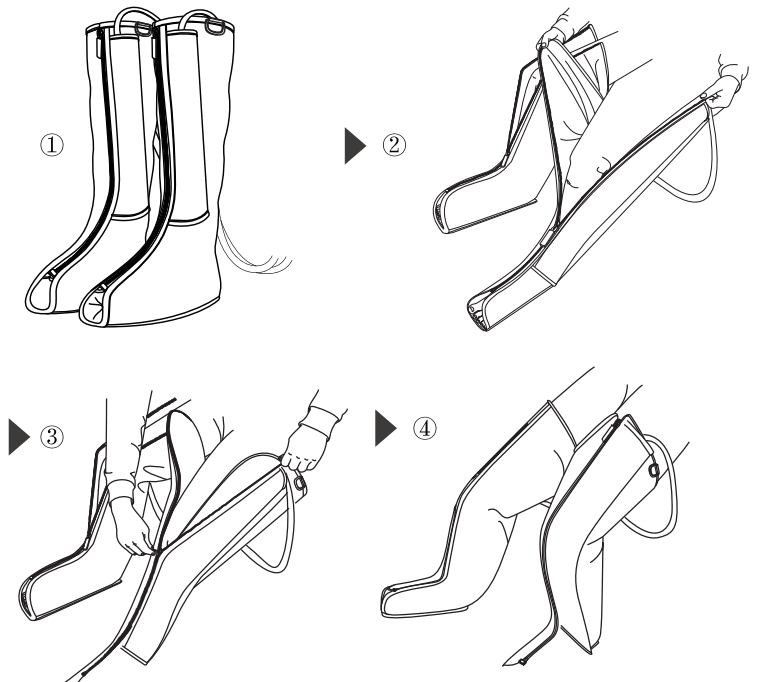Specifications

- Adapter: AC Input: 100~240 Volt , AC 50/60Hz, DC Output: 12V 2A
- Rated power: 24 watt
- Pressure range: 40 - 200 mmhg
- Treatment time: 15 - 60 minutes
- Battery power: Up to 4 hours
- Sleeve size: Inseam: 86.1cm

FAQ


Q1: How dose this product work?

A1: This is an air compression device, 4 air chambers(foot, lower calf, upper calf and thigh) that inflates and deflates with a sophisticated sequence to stimulate hand kneading and stroking, maximizes blood flow throughout the body to help with reduced muscle soreness and/or pains to the treated areas, and improve circulation.


Q2: How many massage modes, and what's the difference?

A2: There is 3 massage modes, please check below for further information:


Mode1: Sequence Mode

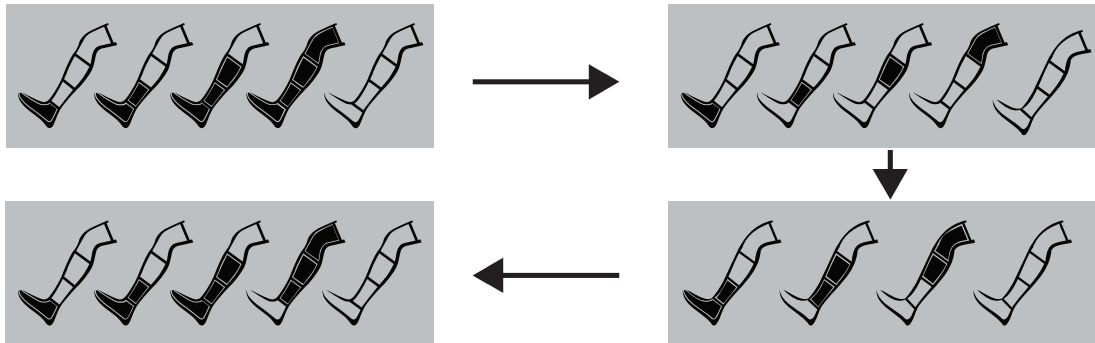Mode 2: Circulation Mode

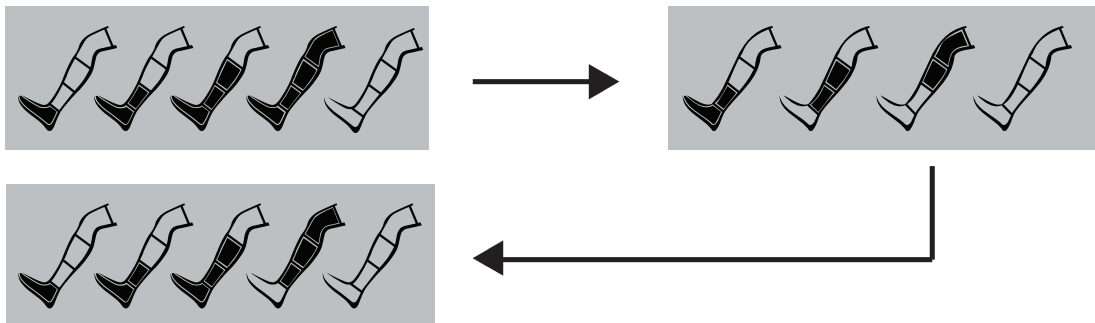Mode 3: Combination Mode

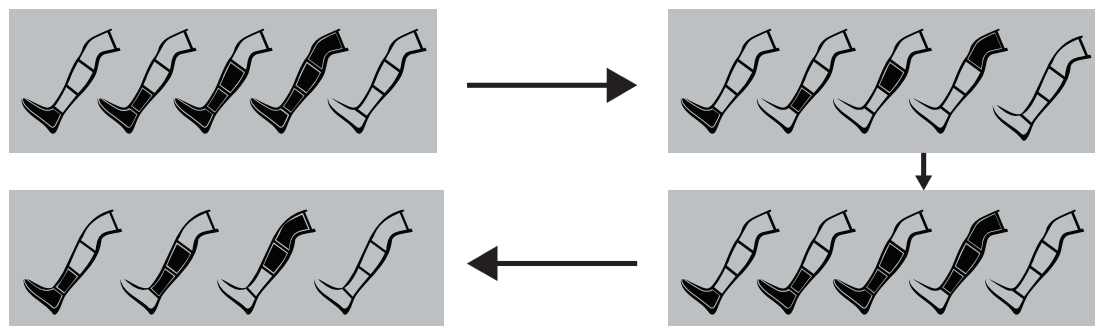Q3: What should I do if I feel the massage strength is too week or too strong?

A3: The intensity is from 40mmhg to 200mmhg, please select the intensity which fits you. If you still feel it's too strong at the lowest setting, please stop using the device, we suggest to consult your doctor before further use.


Q4: How long should I use it?

A4: It is recommended to use it once or twice a day, 15 - 30 minutes at a time. But you can use it longer if your body allows and you don't feel any discomfort. Listen to your body.


Main unit

Find a comfortable position sitting, reclined, or lying down. Press the power button  to start.


- 1). There is a 'Di' sound, the screen will be light up and display current setting.
- 2). It will start with the default setting of: Mode 1 (Sequence Mode), intensity 60mmhg, 30 minutes.

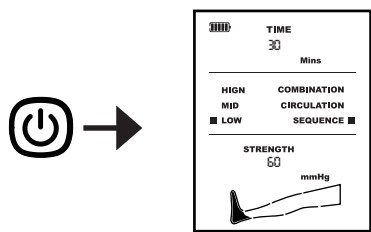Press mode button  to change the message mode. 3 different modes. Please see FAQS A2 for further information.

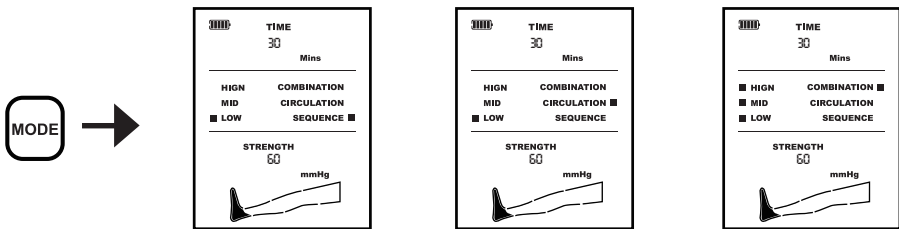
Press zone boost  if you want to fill one specific chamber. Notice: no other chambers will fill.

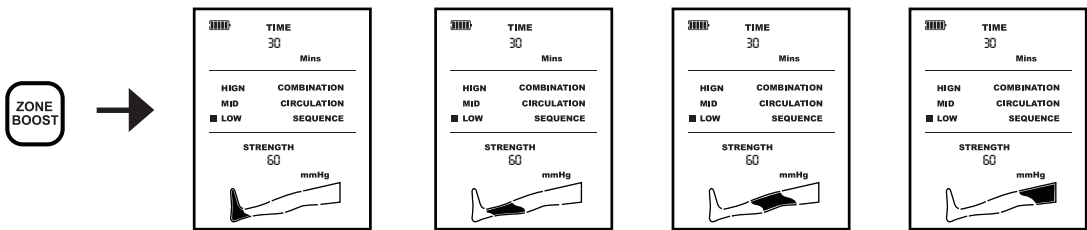Press intensity buttons   to adjust between 40mmhg - 200mmhg.

*We suggest you start with low intensity for the first use, gradually increase after you get used to it.

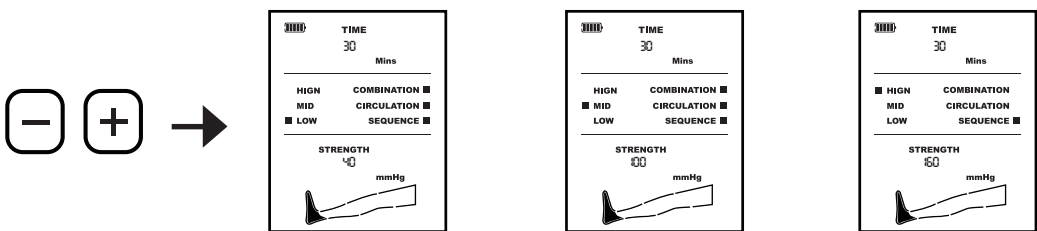IG @Recoverym8 - fb.com/Recoverym8

RECOVERYM8

60 61 23 58

CVR: 41239379

WWW.RECOVERYM8.DK